01306

DE

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ab Seite 3

EN

Instruction manual

starting on page 17

FR

Mode d'emploi

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NL

Handleiding

vanaf pagina 45

Z 01306 M DS V1 0217

Inhalt

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Content

Explanation of Symbols _____	17
Intended use _____	18
Safety Instructions _____	19
Items Supplied _____	22
Device Overview _____	23
Before Initial Use _____	24
Changing the Roller _____	25
Use _____	26
Cleaning and Storage _____	28
Troubleshooting _____	29
Technical Data _____	29
Disposal _____	30

Explanation of Symbols

Safety instructions: Read these carefully and follow them to prevent personal injury and damage to property.

Supplementary information

Circuit symbol for direct current

Read operating instructions before use!

Dear Customer,

We are delighted that you have chosen this **callus remover**.

Before using the device for the first time, please read the instructions carefully and keep them for future reference. If the device is given to someone else, these instructions should also be handed over with it.

The manufacturer and importer do not accept any liability if the information in these instructions is not complied with.

If you have any questions about the device and about spare parts/accessories, contact the customer service department via our website:

www.service-shopping.de

Intended use

- This device is intended to be used to remove calluses as well as dry and wet skin from your feet.
- Only use the device on your feet; do not use it on other parts of your body.
- The device is designed for personal use only and is not intended for commercial applications.
- The device should only be used as described in these instructions. Any other use is deemed to be improper.
- The following are excluded from the warranty: all defects caused by improper handling, damage or unauthorised attempts at repair. This also applies to normal wear and tear.

Safety Instructions

Health-Specific Instructions

- **If you are unsure whether the device may be used, consult a doctor.**
- Do not use the device on injured, infected, loose, or (sun-) burnt skin or in the direct vicinity of acne, birthmarks, warts or varicose veins.
- The device is not suitable for people who suffer from haemophilia, diabetes or an immunodeficiency disorder.
- People who have a pacemaker or belong to another risk group should consult a doctor before using the device.
- If you have skin problems on your feet, consult a doctor before using the device.
- For reasons of hygiene, a roller should only be used by one person.
- If you experience an unpleasant sensation while you are using the device, stop applying it immediately!
- Consult a doctor if skin irritations occur.

Danger of Injury

- **Danger of suffocation!** Keep children away from the packaging material.
- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and/or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.

- Children must be supervised to ensure that they do not play with the device.
- Do not use the device if it is damaged. Before each use, check whether the device and the roller are in good condition. If there is any damage or wear or sharp edges, your skin may be injured.
- Do not use the device on tender, soft skin.
- From time to time, check the area you are treating to avoid penetrating your skin too deeply and causing injuries. Do not hold the roller on one spot for longer than 2 to 3 seconds.
- Make sure that the rotating parts of the device only come into contact with the calluses to be treated. Never touch the roller with your fingers!
- During use, keep the device away from head hair, eyelashes, eyebrows as well as from clothing, brushes, wires, cables or shoelaces to avoid any risks of injury.

Danger from Batteries

- Batteries can be deadly if swallowed. You should therefore keep batteries and the device out of the reach of small children. If a battery is swallowed, seek medical assistance immediately.
- Before inserting the batteries, check whether the contacts in the device and on the batteries are clean, and clean them if necessary.
- Use only the battery type specified in the technical data (see the “Technical Data” chapter).
- Always replace all batteries together. Never use batteries of different types or brands, or batteries with different ca-

capacities. When replacing batteries, ensure the correct polarity (+/-).

- Take the batteries out of the device when they are dead or if you will not be using the device for a prolonged period. This will prevent any damage that might occur if a battery leaks.
- If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you touch battery acid, rinse the affected sites with abundant fresh water immediately and contact a doctor as soon as possible.
- Remove a dead battery from the device immediately. Clean the contacts before new batteries are inserted.
- The batteries must not be recharged or reactivated by other means, dismantled, thrown into fire, immersed in fluids, or short-circuited.

Avoiding Damage to Material, Property and the Device


- **CAUTION! The device may become wet, but must not be fully immersed under water!**
- Use the device only with the roller inserted correctly!
- Only ever use the rollers from the manufacturer. Use the supplied rollers only in this device.
- Do not insert any objects into the openings of the device.
- Protect the device from heat, naked flames, direct sunlight, extreme temperatures, persistent moisture and collisions.
- Do not allow the device to operate unattended. Switch off the device and put on the protective cap when you are not using it.
- In the event of damage/faults, switch the device off immediately.

- Always check the device for damage before putting it into operation. If the device or the rollers display visible signs of damage, the device must not be used. Do not try to repair the device yourself or to replace any individual parts. In the event of any damage, please contact our customer service department.

Items Supplied

- 1 x basic device with protective cap
- 1 x roller set (2x coarse and 2x fine)
- 2x 1.5 V battery, AA type

Device Overview

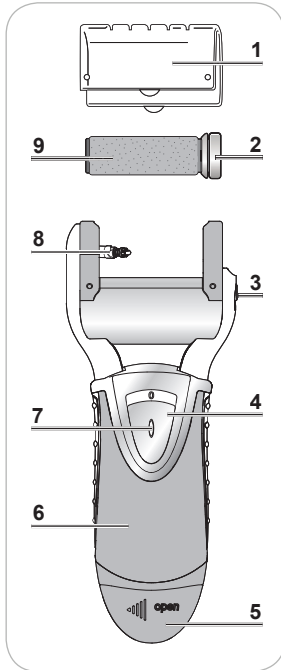- 1 Protective cap
- 2 Roller cap
- 3 Release switch for replacing the roller
- 4 On/Off switch
- 5 Battery compartment cover
- 6 Grip
- 7 Release button of the On/Off switch
- 8 Spindle (for receiving the roller)
- 9 Roller

Before Initial Use

CAUTION!


- **Danger of suffocation!** Keep children and animals away from the packaging material.
-
- Unpack all parts and check the items supplied for completeness (see the “Items Supplied” chapter) and possible transport damage. Remove any stickers and films. If you find any damage to the components, do not use the device (!), but contact our customer service department.
 - Insert the batteries supplied (see the “Inserting / Changing Batteries” section).

Inserting / Changing Batteries

CAUTION!

- Replace the batteries when the power of the device starts to wane.
- Switch off the device before you replace the batteries (move On/Off switch to position **0**).

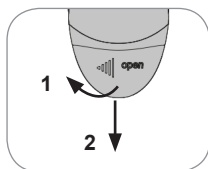1. Rotate the battery compartment cover a quarter turn in the direction of the arrow.
 2. Then take off the battery compartment cover.
3. If necessary, remove the used batteries and dispose of them in an environmentally friendly way.

4. Insert two new 1.5 V batteries of the AA type. Pay attention to the correct polarity (+/-).
5. Place the battery compartment cover back on the battery compartment and turn it counter to the direction of the arrow to shut it.

Changing the Roller

CAUTION!

- A damaged/worn roller can cause injuries to the skin and must be replaced.
- Make sure that the device is switched off before changing the roller.

Depending on the thickness of the callus, the appropriate roller (9) must be used. To remove thick calluses (e.g. on your heels), use the coarse black roller. To remove dry and thin callus areas, use the fine grey roller.

1. Remove the protective cap (1).

2. Hold the device so that the roller cap (2) points upwards and slide the release switch (3) downwards. The roller (9) will "spring" out.
3. Pull the roller (9) out of the holder.
4. Using your finger, press the new roller (9) out of the protective sleeve and hold the roller cap (2) between the fingers of your other hand.

5. Slide the narrow end of the roller (9) into the holder in the device until you hear it engage.

6. Check that the roller (9) fits securely by shaking the device gently from side to side.

Use

CAUTION!

- Check the device/roller for damage and wear before each use.
- Make sure that the release switch for replacing the roller is not accidentally pressed during use!
- Only use the device on uninjured calluses on your feet.
- If as you are using the device you experience any pain or skin irritations, stop using the device immediately!
- Do not allow the device to rotate on the same spot on the skin for longer than 3 seconds. From time to time, check the area you are treating to avoid penetrating your skin too deeply and causing injuries.
- Do not press the roller against the skin with force. The skin could be injured.
- To prevent accidents, keep the device away from head hair, eyelashes, eyebrows, textiles, cables and shoe laces etc.

- For dry application, make sure that there are no cream or care oil residues on your feet.
- The device may be used under the shower and in the bathtub. **CAUTION! The device may become wet, but must not be fully immersed under water!**

Removal of Calluses

1. Remove the protective cap (1) from the device.
2. Switch on the device by simultaneously pressing the release button (7) on the On/Off switch (4) and sliding the On/Off switch (4) upwards (I will now be visible beneath the On/Off switch (4)).
3. Place the rotating roller (9) against the callus which is to be removed and move the device back and forth and to the sides on the callus. Move the device slowly and without pressure.

- If the device is pressed too firmly against the skin, the roller (9) will stop automatically to prevent skin injuries.
4. Switch off the device after use by simultaneously pressing the release button (7) of the On/Off switch (4) and sliding the On/Off switch (4) downwards (0 will now be visible above the On/Off switch (4)).

After Treatment

1. After treatment, apply cream to the feet so that the skin becomes soft and tender.
2. Clean the device after each use (see the “Cleaning and Storage” chapter).

Cleaning and Storage

CAUTION!

- Make sure that the device is switched off before cleaning it!
- **The device may be rinsed under running water, but must not be fully immersed under water! It could be damaged!**
- Only use water to clean the rollers.
- Do not use any corrosive or abrasive cleaning agents to clean the device. These can damage the surface.

1. Take the roller (9) out of the device (see the “Changing the Roller” chapter) and brush off any flakes of skin.
2. Hold the head of the device under running water for 5 to 8 seconds. Also clean the roller (9) under running water.
3. Allow the device to dry thoroughly and allow the roller (9) to dry completely in the air.
4. Insert the roller (9) back into the device and put on the protective cap (1).
5. Store the device and accessories in a dry place which is not accessible to children and animals.

Troubleshooting

If the device does not work properly, first check whether you are able to rectify the problem yourself.

Do not attempt to repair a defective electrical device yourself!

Problem	Possible cause / solution
The device is not working.	Are the batteries too weak or incorrectly inserted? If necessary, insert new batteries or check the polarity (see the “Inserting / Changing Batteries” section).
The device stops during treatment.	Is too much pressure being exerted? For safety reasons, the device stops automatically when too much pressure is applied on the skin.
The calluses are not being properly removed.	The roller (9) is worn. Replace the roller (9) (see the “Changing the Roller” chapter).

Technical Data

Article number: 01306
Model: AP-3WD
Voltage supply: 3 V DC (2 x 1.5 V battery, AA type)

Disposal

The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner and take it to a recycling centre.

Dispose of the device in an environmentally friendly manner at a recycling centre for used electrical and electronic devices. It should not be disposed of in the regular household waste. For more information, contact your local authorities.

The batteries should be removed prior to disposal of the device and disposed of separately from the device. To protect the environment, non-rechargeable and rechargeable batteries may not be disposed of with normal household waste, but must be taken to suitable collection points. Also note the relevant statutory regulations governing the disposal of batteries.

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(Calls to German landlines are subject to charges. The cost varies depending on the service provider.)

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All rights reserved.
DS Produkte GmbH, Am Heisterbusch 1