CIRCULATIONS- TRAINER

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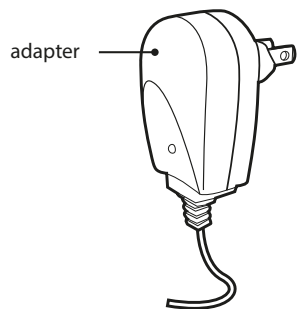
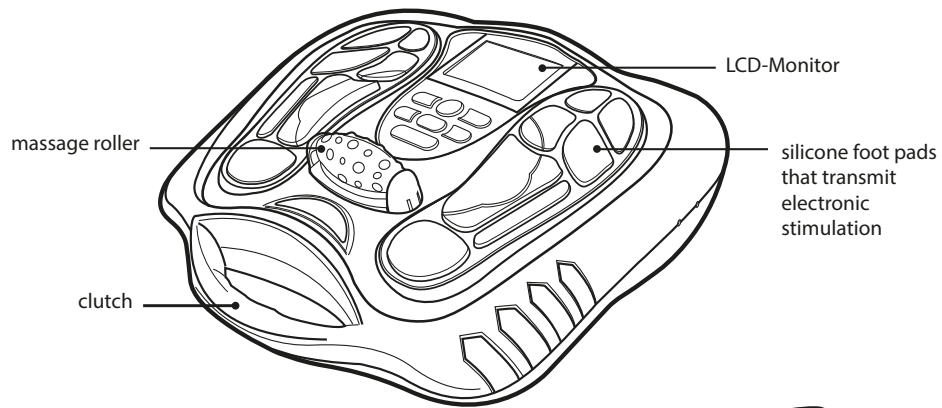
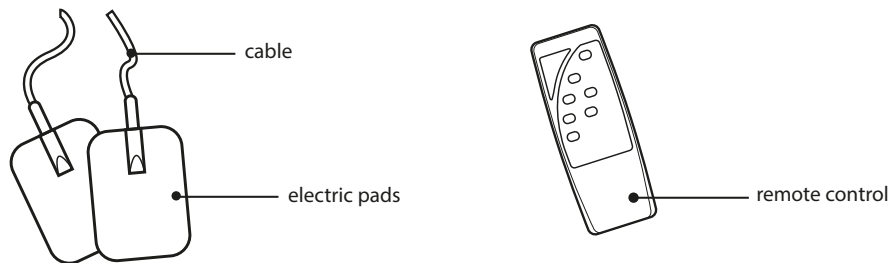
Dear customer, dear customer,

we are pleased that you have chosen this product. If you have any questions, contact customer service via our email address:

Kundenservice@hsphanseshopping.de

i Please read the manual carefully before using the device for the first time and keep it in a safe place. These instructions must be handed over with the device if it is passed on. The manufacturer and importer assume no liability if the information in this manual is not observed! In the context of constant further development, we reserve the right to change the product, packaging or accompanying documents at any time.

 All safety instructions are marked with this symbol. Read these carefully and follow the safety instructions to avoid personal injury and damage to property to avoid.

 All warnings are marked with this symbol. Read these carefully and follow the warnings to avoid personal injury and damage to property to avoid.

 Disposal: Do not dispose of electrical appliances as unsorted.

 The device may only be used in dry rooms.

Summary Application

Operation Circulation Maxx:

1. Remove the Circulation Maxx and accessories from their packaging.
2. Place batteries in the remote control. The unit also works without remote.
3. Connect the device with the supplied cable.
4. Switch on the Circulation Maxx. Press the „On/Off“ button for three seconds.
5. The LCD monitor should light up blue.
6. Remove your shoes, socks or pantyhose and place your bare feet on the intended places on the device. Both feet must be on the pads for the Circulation Maxx to work.
7. PLEASE NOTE THAT THE AIM IS NOT TO REACH LEVEL 99! YOU WILL SIMPLY FIND A COMFORTABLE LEVEL FOR YOU.

FAQ

Q: I cannot feel any vibration of the Circulation Maxx.

A: The device does not vibrate. If you have placed both feet on the pads and selected the appropriate intensity level, you should feel the electrical stimulation.

Q: I cannot feel any change in my feet during use.

A: You may have chosen too low an intensity level. Adjust it using the arrow keys.

Q: Yesterday I used the device at level 45, today I can't feel anything at this level. A: Your body's response to the application can change daily. On some days you should choose a higher intensity level to optimize the therapy.

Q: My legs hurt after use.

A: Either you have used Circulation Maxx too high or too often in the last few days. Give your legs time to relax and then continue with the application.

Q: I feel no stimulation in my feet even at a high intensity level.

A: You may be dehydrated. You should drink enough to support the electrical impulses.

How does electronic stimulation work?

field of application: Medical treatment purpose

Please note that it can take up to four weeks for you to notice initial results. We recommend using Circulation Maxx at least once a day. Finish the application of the device, you will also decrease the blood circulation effect feel it. Remember to drink plenty of water because that will increase the benefits.

The electric pads are intended for the treatment of hard-to-reach areas. You can use the electric pads and the Circulation Maxx at the same time.

Electronic nerve stimulation works like a massage that relieves pain, promotes circulation, relaxes stiff muscles, swollen feet and ankles. The nerves are stimulated by the electronic stimulation in the feet and the places where you attach the electric pads.

Suitable for: Please read the safety information before using the device.

Electronic nerve stimulation is an unobtrusive, safe form of nerve stimulation designed to relieve pain. The Circulation Maxx uses proven neuromuscular electronic stimulation therapy that sends regular micropulses through the soles of the feet. This type of electrical nerve stimulation has been clinically proven to be safe, effective and can be applied in the comfort of your own home.

Operation

1. Place your bare feet on the intended places on the Circulation Maxx (left foot on the left pad and right foot on the right pad).
2. Turn the power on by holding down the On/Off button for three seconds. The LCD monitor will turn blue.
3. Select the desired mode by using the „Mode“ ▲▼ arrow keys. You can choose from 25 different modes.
4. Increase or decrease the intensity with the „Sole“ ▲▼ arrow keys (stimulation of the soles of the feet) or the „Body“ keys (stimulation by the glued-on electric pads).
5. The intensity and modes can be changed at any time. However, if you switch off the device, the settings will be reset.
6. To end the treatment, simply press the „On/Off“ button again for three seconds.

FAQ Electro-Pads

- **No stimulation** Are the cables connected correctly? Have you removed the protection from the electric pads?
- **Low stimulation** Do the electric pads stick correctly to the skin? Are the electric pads dirty? Are the electric pads damaged?
- **The skin reddens and hurts** Do you use the pads for too long? Are the electric pads too dry? Do the electric pads stick too strongly to the skin? Are the electric pads dirty? Are the electric pads damaged?
- **The stimulation is interrupted** Have the electric pads become detached from the skin? Has the cable come loose from the pads? Is the treatment time over?

Safety instructions for electric pads

- Never glue two electric pads together.
- Make sure that the electric pads are always clean and keep them away from direct sunlight and high temperatures.
- If the electric pads do not stick sufficiently or are dirty, clean them with a damp cloth or use new ones. Never clean the electric pads with chemicals.
- To protect the electric pads, always stick them back onto the pad protectors after treatment.

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Clean the device after use. Do not immerse the device in liquids. Do not use petrol, window cleaner, glass cleaner, furniture polish or paint thinner. Do not use a brush or fingernails to clean the electrodes to avoid damaging the electrode surface. Store the electrode gel pads on the supplied plastic pad protectors.

Safety instructions

Persons with the following characteristics should consult a physician before using Circulation Maxx.

1. acute diseases
2. malignant tumours
3. infectious diseases
4. pregnancy
5. heart failure or disorder
6. high fever
7. abnormal blood pressure
8. skin problems
9. in medical treatment

Attention

If the device does not function correctly or you notice irregularities, please stop the treatment.

If you want to move the electric pads to another part of your body, make sure to disconnect them from the power supply first.

Do not attempt to attach the electric pads to another person during treatment. Severe electric shocks could result.

Do not use the device with children or persons who have your own will. I can't express. Do not use the Circulation Maxx if you are younger than 18 years.

Do not use in rooms with high humidity, such as bathrooms. Do not use the product while taking a bath or shower. You may experience a severe electrical shock that may result in serious injury.

Do not use the device while you sleep. Irregularities may occur with the Circulation Maxx or the electric pads may slip, resulting in health hazards.

Do not use the device while driving. Electronic nerve stimulation could distract you and cause an accident.

Remove the electric pads from the skin after use. Skin irritation may occur if the electric pads are left to adhere to the skin for too long.

Make sure that the electric pads do not come into contact with metallic objects such as belt buckles or chains.

Do not use a mobile phone or other electronic device right next to the device on.

Do not use the device near the heart, around the neck, on the head, around the neck, on the neck, on the head, on the neck, on the head or on diseased skin. In these places it may be necessary during application to accidents or health hazards.

If you place the electric pads close to the heart, it may be dangerous, especially if you have Susceptibility, lead to irregularities.

Do not use Circulation Maxx at the same time as other medical devices or directly on ointments.


Please use the Circulation Maxx only for the purposes described in this manual.

Please use the cables and plugs only in the intended locations. Otherwise accidents or electric shocks may occur.

Please do not dismantle the device or carry out any modifications. This could cause fire, short circuits or accidents.

Description of the different modes
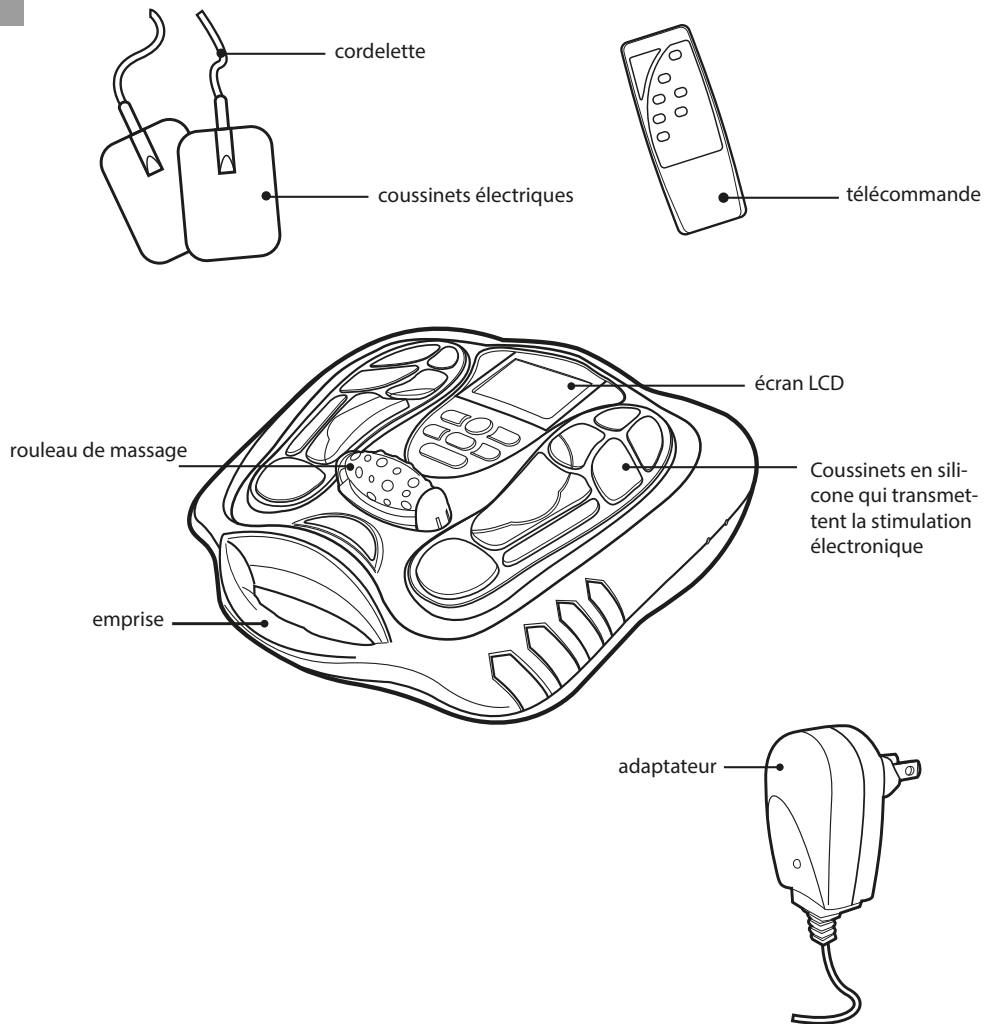
| Mode | Effect | Mode | Effect |
|------|----------------------------------|------|-------------------------|
| 1 | Pressing the reflex zones | 14 | Scratching and pressing |
| 2 | reflex zones | 15 | Deep drawing |
| 3 | Kneading the reflex zones | 16 | Shaping |
| 4 | Knocking the reflex zones | 17 | Easy hitting |
| 5 | scratches | 18 | Massage and tapping |
| 6 | presses | 19 | Press |
| 7 | Massage | 20 | Reaming and rolling |
| 8 | Pressing massage | 21 | Presses |
| 9 | Pressing and pressing | 22 | light impacts |
| 10 | Pressing the reflex zones | 23 | Reflex zone massage |
| 11 | Easy beating of the reflex zones | 24 | Shiatsu |
| 12 | Kneading | 25 | Rolls and kneading |
| 13 | drums | | |

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HSP Hanse Shopping GmbH